




It is our job to make
your glasses
comfortable

Hasegawa Bicoh co., ltd.

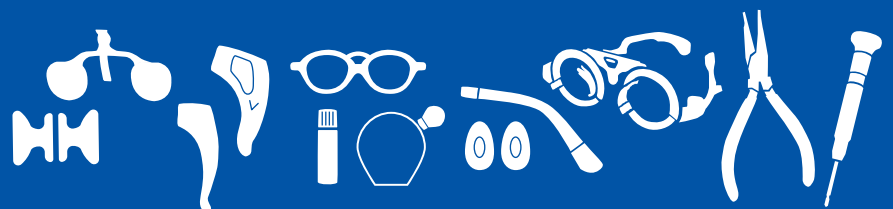

CONTENTS

1 CELL SEAL U
CELL SEAL M

2 MEGALOCK V
MEGALOCK F
MEGALOCK M

3 MEGALOCK BALANCE

4 MEGANE GRIP
UP & DOWN NOSE PAD

5 ROUGLASS L

6 CLEAN COAT

7 BiBa

8 SILICONE NOSE PAD DX
CUSHION NOSE PAD C-one
CUSHION HILON NOSE PAD

9 HILON NOSE PAD EX
HILON NOSE PAD S
HILON NOSE PAD R
CUSHION NOSE PAD AR

10 KUBI FREE

11 DRIVER(+) FOR REMOVING
BROKEN SCREWS

12 GACHIDORA

13 GACHIDORA ZARA

14 CSL PLIER

15 CS PLIER
CELL MORI
SILICONE COAT SET

16 SILICONE TEMPLE TIP (ROUND HOLE)
SILICONE TEMPLE TIP (FLAT HOLE)
SILICONE TEMPLE TIP (SLIM TYPE)

17 SOFURE TEMPLE TIP
SILICONE CURL TIP FOR ADULT
SILICONE CURL TIP FOR KIDS

UNIVERSAL TRIAL FRAME

18 580UniV Pro
580UniV Pro – LO

20 MC674
MC486

FIX TYPE TRIAL FRAME

21 SIMPLE BC
SIMPLE BC – LO
COLOR SIMPLE BC

22 FITTING MEASURE
TRIAL SET FOR HIGH CURVE

23 CLUTCH GLASS α
for “EYELID PTOSIS”


All made in Japan



To prevent slip down of the plastic frame

CELL SEAL® U

Point for improvement

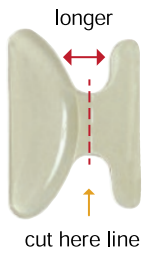

- Just paste on the nose pad.
- When the nose pad is raised, the distance between the eyes and the lens is separated. Eyelashes and cheeks are less likely to touch the eyeglasses.
- The shape of the cell seal has been changed to make it easier to turn up. Easy to peel off and paste by making cuts in the release paper.
- It is transparent color, made of silicone which is safe to the skin.



<PAT no. 3157348>

How to use

Raise the nose pad



Paste the nose pad of the frame sandwiching.

Narrow the nose width




Cut and paste the reinforcement part.

Nose pad part Cut Reinforcement part

Size (Thickness)	S (1.3mm) / M (1.8mm) / L (2.5mm) / LL (3mm)			Package	1 pair
Material	Silicone			Box	Including 10 pairs.

Small Size CELL SEAL

CELL SEAL® M

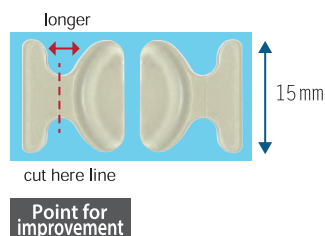

<PAT no. 3157348>

- For small plastic frame.
- Easy to turn up and paste.

Size (Thickness)	M (1.8mm) / L (2.5mm) / LL (3mm)			Package	1 pair
Material	Silicone			Box	Including 10 pairs



Suitable for various size temple tip

MEGALOCK® V



Black

Blue

Green

Pink

Actual size

* Width of temples : about 5mm - 10mm

- Just insert to the temple tip. Glasses are fixed.
- Especially recommended for those who like to move such as sports.
- Since the distance between the eye and the lens is fixed, the burden on the eye is reduced.

Include adjust parts



For Glasses with thin temple tip.

Package	1 pair	Include adjust parts & case
Color	BLACK / BLUE / GREEN / PINK	
Material	Silicone	
Box :	1) 1 color : 10 pair 2) Assortment : 10 pair (Black:3 / Blue:3 / Green:2 / Pink:2)	



Large size Megalock

MEGALOCK® F



Black

Brown

Beige

Actual size

* Width of temples : about 4mm - 8mm

Package	1 pair	Include adjust parts & case
Color	BLACK / BROWN / BEIGE	
Box	10 pairs	
Material	Silicone	



Small size Megalock

MEGALOCK® M



Black

Brown

Green

Actual size

* Width of temples : about 3mm - 6mm

Package	1 pair	Include adjust parts & case
Color	BLACK / BROWN / GREEN	
Box	10 pairs	
Material	Silicone	





NEW

Products that reduce the weight of nose.

MEGALOCK® BALANCE

**Attention** BALANCER "Weight" inside.

Black



Beige



Gray



Actual size


* Width of temples : about 5mm - 10mm

- A weight is built into the conventional MEGALOCK.
- When worn towards the tip of the temple, the center of gravity of the glasses moves backwards, reducing the burden on the nose.
- When worn behind the ear, it can also be used to secure glasses in place, so it is recommended for people who engage in physical activity such as sports.

Lock

1 Securing them behind your ears.

Lock your glasses.

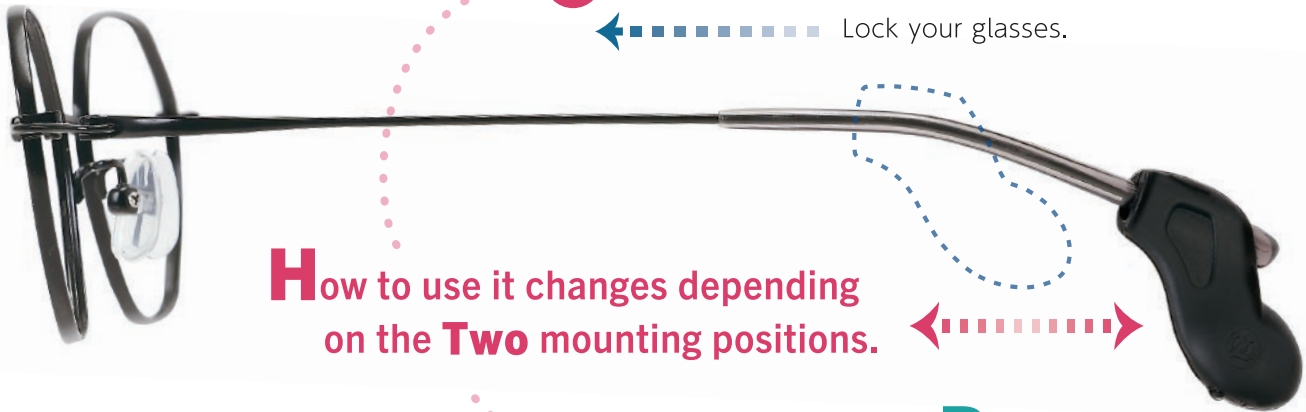How to use it changes depending on the **Two** mounting positions.

Relax

2 Attach to the tip of the temple.

Reduces the burden on the nose by approximately 20% by the weight.*

※ Compared to our company.



<PAT no. 3242259>

Include adjust parts



For Glasses with thin temple tip.

Package	1 pair	Include adjust parts
Color	BLACK / BEIGE / GRAY	
Box	8 pairs	
Material	Silicone, Brass	



Just cover the temple tip by the silicone tube.

MEGANE GRIP

- Just insert to the temple tip.
- Suitable for temple tip because stretch silicone.
- Safety silicone material.




Package	1 pair	Box	10 pairs
Color	Brown	Material	Silicone
Size	M / L	Outer/inner	M : 2.6 ϕ / 4.3 ϕ L : 3.1 ϕ / 4.8 ϕ

For the progressive lens use

UP & DOWN NOSE PAD[®]



- The nose pad moves up and down by 5 mm. It makes it easy to see nearby.
- V character shape. Vertically up and down. The nose width does not change.
- Silicone surface

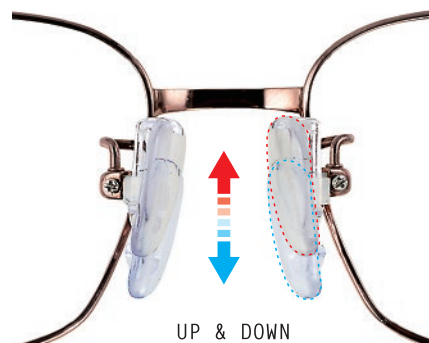

<PAT no. 4833954>

		Display set	
Package	1 pair	5 pair set	12 pair set
Type	Screw in / Push on	Screw in	Screw :10 / Push :2
Size	15.5 × 8mm	Hole size	1.2 ϕ
		Material	Nylon, Polyethylene


Since the near portion of the lens approaches the eye, You can see small characters and others nearby.



Change to bifocals just by attaching to glasses.

ROUGLASS® L

Just pinch the clip and insert it inside the eyeglasses. Add the lens power to make it easy to see nearby.

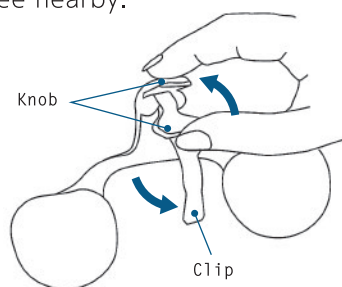

- It can be installed with one touch.



Sales point

- I can not see detailed characters clearly.
- For nearness a little bit hard to see, but a bit more.
- I want to use my current glasses.
- It is troublesome to carry reading glasses on a trip or going out.


Reference for lens power

AGE	POWER
40-50	+1.00
45-55	+1.50
50-60	+2.00
55-65	+2.50
60-70	+3.00
70~	+3.50



<The licence of "manufacture and sales of medical equipment" no. 13B3X00177>

Power	+1.00 / +1.50 / +2.00 / +2.50 / +3.00 / +3.50
Material	Lens : Acrylic Clip : Nylon
Display	9 pieces(+1.00 :1 / +1.50 :2 / +2.00 :2 / +2.50 :2 / +3.00 :1 / +3.50 :1)

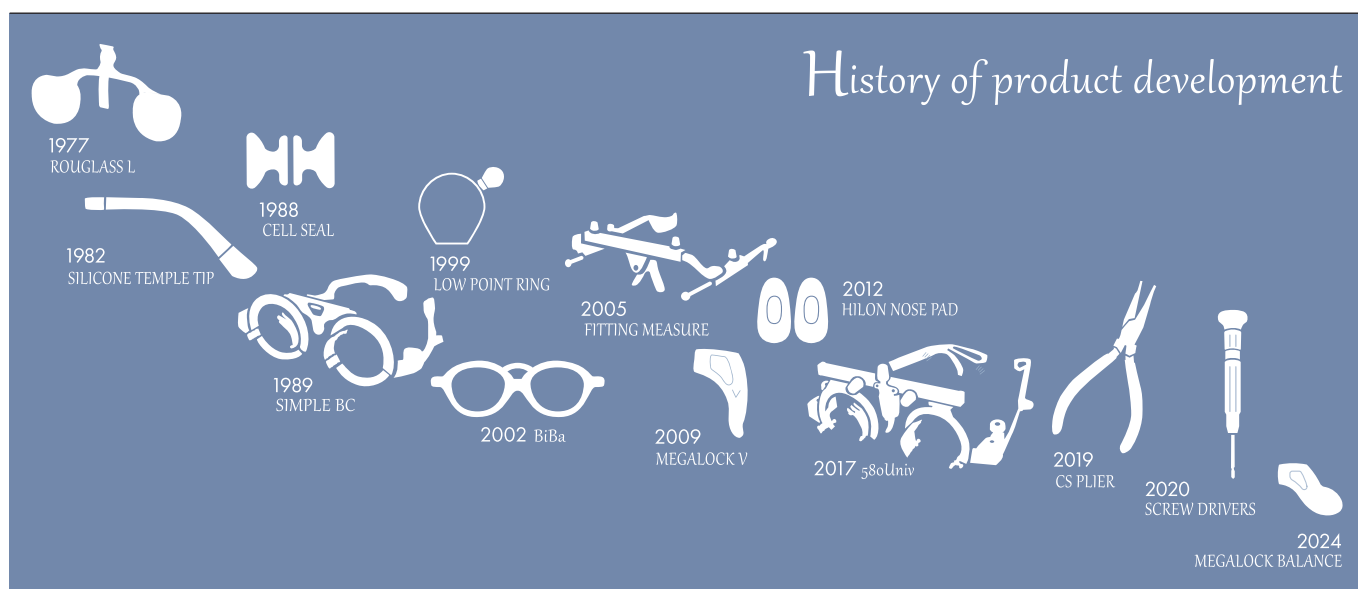



Fluorine blending in lenses, slippery, long-lasting glasses

CLEAN COAT

Discontinued

Make a coating film of fluorine and protect glasses from scratches and dirt.

Also for smart phones, PC monitors and other uses other than glasses.



Bulk	10mℓ
Material	Surfactant · Fluorine oil

Both cleaners and coatings can be done!
Alcohol is not used at all.

Effect on lens

- Removes a wide variety of dirt on the entire lens surface.
- By making a coating film of fluorine on the whole lens surface
 - Prevent contamination and scratches.
 - Easily remove dirt.
 - Touch feeling when wiping the lens refreshing It is slippery.

Effect on frame

- By making a coating film of fluorine on the whole frame
 - It protects from hair tonic for Temple and temple end.
 - If you are allergic to the rash in Temple Prevent.
 - Prevention of whitening of plastic and tortoise shells.









Package	For sale : 10mℓ (20 pcs.)	For store demonstration : 200mℓ (3 pcs.)	Replenishment : 1ℓ (1 pcs.)
---------	------------------------------	---	--------------------------------



Optical frame for a baby (1-2 age)

BiBa









Emphasis on safety · cuteness, hold tolerance for young children.
 Since after parts are also available, various customization is possible.
 Please try the shape that you can feel comfortably without disliking glasses.

	
①+②+③+⑦	①+②+④ / ⑤+⑥
Head band type	Curl tip type
	
①+②+④ / ⑤+⑥+⑦	①+②+④ / ⑤+⑦×2本
Curl tip & Head band type	T-style type
	
①+②+④ / ⑤+⑦×2本	①+②+④ / ⑤+⑦×3本
Goggle type	X-style type



Size	Frame : 36 / 38
Color	Pink / Blue

After part

① Front frame	② Silicone nose pad	③ temple for the band	Parts set
			
④ short temple ⑤ Long temple	⑥ Curl tip	⑦ Band	
		Front  Back 	

Color : Pink / Blue



Standard of the silicone nose pad

SILICONE NOSE PAD DX




Size	L (18mm)/M (16mm)/S (13mm)		
Color	Pad : Transparent Metal : Gold or Silver	Hole size	1.2 φ
Package	1Box : 9 pairs (Gold :6 pairs/Silver :3 pairs) (each size) Bulk : 50 pairs (each size)		

Gently absorb the movement of eyeglasses by the cushion function

CUSHION NOSE PAD C-one

- A resin cushion mechanism is provided in the conventional core portion. By integrally molding (making all shapes with one material) It enables unprecedented strength, springiness and flexibility We realized comfortable wearing feeling.
- Newly developed material Using ethylene elastomer, it is excellent in durability and chemical resistance, It is also used for medical use (clear food sanitation law), it is also a gentle pad on the skin.

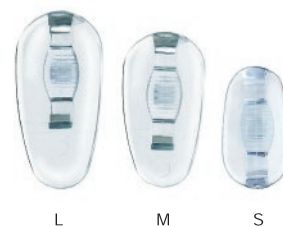


<PAT. no.3208475>

Size	L (17mm)/M (14.6mm)/S (11.6mm)		
Color	Transparent	hole size	1.1 φ
Type	Screw in / Push on		
Package	1 Box :25 pairs (sizes and types each) Bulk :50 Pairs (sizes and types each)		

Clear cushion pad of special plastic

CUSHION HILON NOSE PAD

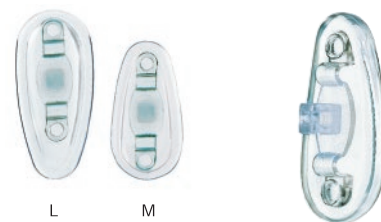


Size	L (17mm)/M (14.4mm)		
Color	Transparent	Hole size	1.1 φ
Package	1 Box :9 pairs (each size) Bulk :50 pairs (each size)		

Hilon nose pad SERIES

- High performance nose pad using special engineering plastic.
- Nose pad traces are difficult to attach to the nose.
- Cosmetics and sebum dirt are difficult to adhere.

There is cushioning property by adoption of silicon ring

HILON NOSE PAD EX

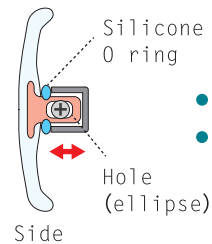

- Soft cushioning.
- No sound and no clattering.

<PAT. no. 3169961>

Size	M (14.4mm)/S (13mm)
Hole size	1.2×1.6mm (1.2 screw available)
Color	Pad :Clear beige / Metal :Gold
Package	1 box :9 pairs (each size) Bulk :50 pairs (each size)



Standard Hilon pad

HILON NOSE PAD S



Size	14.4mm	Hole size	1.2 φ
Color	Pad :Clear beige / Metal :Gold		
Package	1 box :9 pairs / Bulk :50 pairs		

Hilon pad for Men

HILON NOSE PAD R



Size	15mm	Hole size	1.2 φ
Color	Pad :Clear blue / Metal :Silver		
Package	1 box :9 pairs / Bulk :50 pairs		

Cushion function of the new mechanism realized

CUSHION NOSE PAD AR



- Adopted comfortable material "ethylene elastomer".



← This gap is cushion

<PAT no. 4833790>

View from top side.

Size	14mm		
Color	Transparent	Hole size	1.1 ϕ
Package	1 box : 9 pairs / Bulk : 50 pairs		





Pad that moves 360 degrees freely

KUBI FREE

NEW



- There is a “constriction” between the pad and the screw hole, which makes the pad body movable, so it naturally fits any shape of the nose.
- It is good fit for your nose which makes the pad less noticeable.
- There is no material inlet gate in the main body of pad, so it is safe and smooth to your skin.
- Made of elastomer, which is soft but don't crack.




View from the side.



Attention

“Constriction”

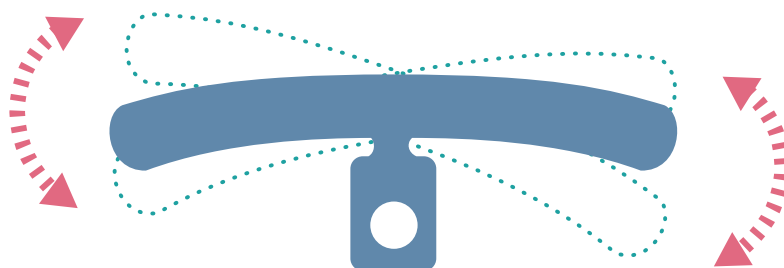

The pad moves due to the effect of the “constriction”, It fits your nose naturally.



<PAT no. 3234977>

Size	14.4mm		
Color	Transparent	Hole size	1.1 ϕ
Package	1 box : 25 pairs / Bulk : 50 pairs		





Easy to remove broken screws

Driver (+) for removing broken screws


NEJI TORU[®] plus ⊕

Bicoh's
Eyeglass tools


- A new screwdriver that can remove broken screws.
- This driver is dedicated screwdriver to remove with a crushed screw head.
- For hinge screw, rim screw, pad screw. Phillips(+) screw only. (Depends on size and condition.)



• Hinge screw

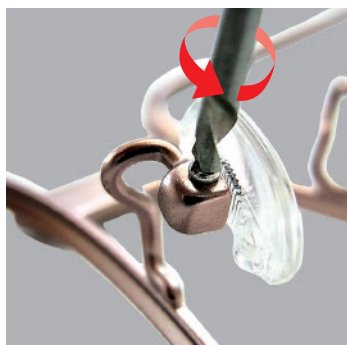

• Pad screw

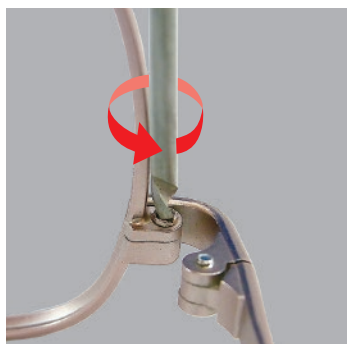

• Rim screw





Actual size

length :122mm
weight:46g

Attention

Try this before sending for repair!

Can remove a crushed screw.

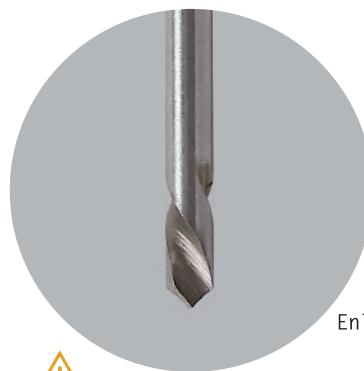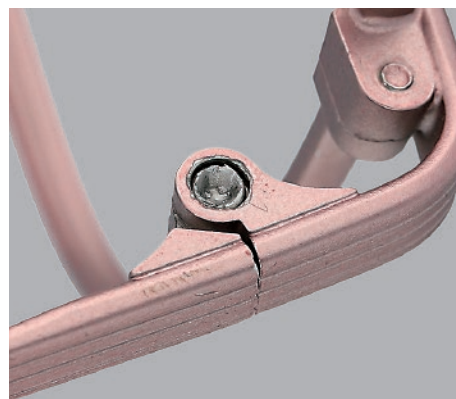

Enlarged image of the tip



This driver is for counterclockwise turn only.
Do not turn clockwise.

Optional replacement parts

bit for the driver



Actual size

Screw diameter	1.0mm 1.2mm 1.4mm
Package	Driver 1 pc. Wrench for replace the driver tip.



Precise phillips(+) driver

GACHIDORA®

Bicoh's
Eyeglass tools



Actual size


length : 115mm
weight : 33g

- A screwdriver that fits very well in the groove of the screw.
- It doesn't move rattlingly.
- The power to turn is efficiently transmitted to the screw.

Attention

Precision Phillips(+) screwdriver


The tip of the screwdriver was designed only for eyeglass screws.

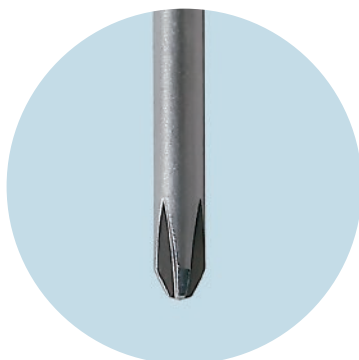


Enlarged image of the tip



• Hinge screw

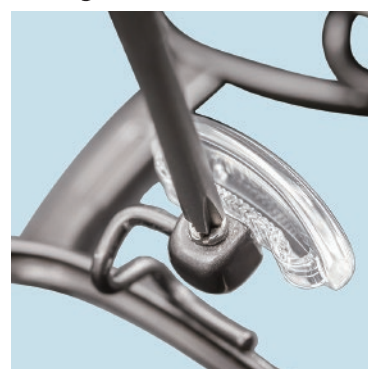

• Pad screw

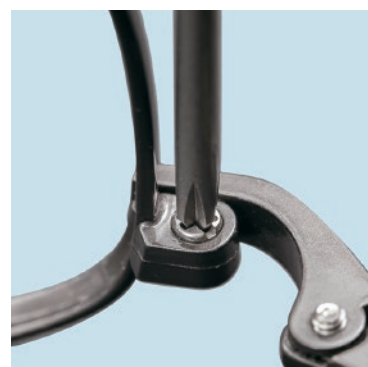

• Rim screw

<PAT no. 3234907>

Screw diameter	1.0mm 1.2mm 1.4mm
Package	Driver 1 pc. Wrench for replace the driver tip.

Optional replacement parts

bit for the driver



Actual size

accessories

Silicone tube for
tip protection to pick up
screws.



Bicoh's
Eyeglass tools



NEW

Precise flat head(−) driver

GACHIDORA® ZARA

- Gachidora's flat head driver is new released.
- By the special "surface treatment", it fits into the screw groove and stop slipping.
- The special "shape" and rounded tip are safe for glasses and hands if slipped while handling.



• Hinge screw



• Pad screw




• Rim screw

Zara

Attention

Non-slip in the groove of the screw.

The tip is diamond-plated,
so stop slipping and safe.

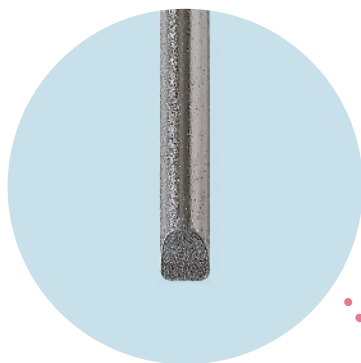

Enlarged image of the tip



Actual size

length :115mm
weight:33g

<PAT no. 3239937>

Screw diameter	1.0mm 1.2mm 1.4mm
Package	Driver 1 pc. Wrench for replace the driver tip.

Optional replacement parts

bit for the driver



Actual size



The plier can adjust various nose pad arms.

CSL PLIER

BicoH's
Eyeglass tools

- CSL plier can adjust "click-in type" nose pad arm smoothly.
- Good hold the metal in the nose pad part with "grooves inside the tool" it can be adjusted with small power.
- It can also be used for "screw-type" pad box.




2020


"The most emotional
eyewear" Grandprix




New Function of CSL plier (CS plier do not have)

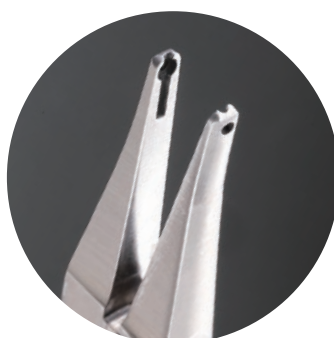


Vertical groove at the tip.

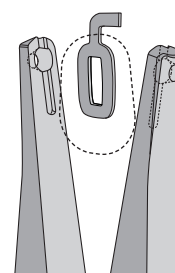

Very useful for adjusting these types of Clings.

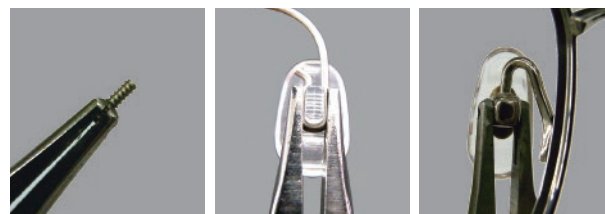

- It also has the function of a CS plier.

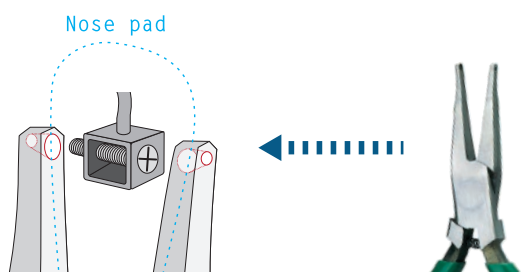

Please grasp the screw so that the screw is in the hole of the plier tip.

<PAT no. 3223800>

Package

Plier 1 pc.
Accessory A (fix for end tip) B (fix for temple)

Accessory parts

Insert to the top of the plier.



A (fix for end piece)




B (fix for temple)





Bicoh's
Eyeglass tools

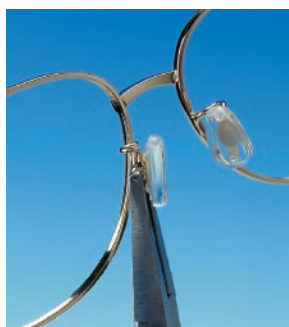

The plier for adjust of nose pad

CS PLIER

- The CS plier is useful tool. Adjust nose pad box easier.
- Since you can grasp only the pad box surely, There is nothing to hurt.
- You can also use it for regular pads. (silicone nose pad, hard nose pad, etc.)
- Include accessory parts. (same as the CSL Plier.)

Package

Plier 1 pc.
Accessory A (fix for end tip) B (fix for temple)

For celluloid
and acetate




No.1

No.2

No.3

Silicone nose pad for the plastic frame

CELL MORI

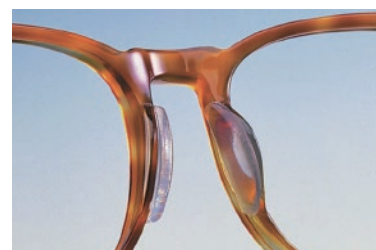

- Please bond with "Silicone Coat Set".
- 1) Take off the oil stains on the part to be coated.
 - 2) When partially painting in a place without a step please mask with a cellophane tape etc.
 - 3) Coat the primer and let it dry for more than 5 minutes.
 - 4) Have the coat part level and apply silicone with a brush.
 - 5) Paste the cell mori, touch lightly with your fingers and pull out the air.
 - 6) Dry about 2-10 hours.

Size	No.1 (1mm) / No.2 (1.5mm) / No.3 (2.2mm)
Package	1 box : 12 pairs (per size)
Material	Silicone



Primer



Silicone glue



Coating brush

For celluloid
and acetate

Functions

- Protect the human skin from plating, apply to the parts to prevent rash.
- Replacement of silicone portion of silicone processed frame, for repair bonding.
- It can also be used for coat for nose pad of plastic frames.

⚠ Please close the lid tightly and store this silicone glue in the refrigerator (it may harden).

Package

Silicone glue×1 pc. / Coating brush×1 pc. / Primer×1 bottle


It can adhesiving and coating

SILICONE COAT SET

- Apply to the nose pad or temple tip part, etc. to prevent pain and gaps.



SILICONE TEMPLE TIP (ROUND HOLE)

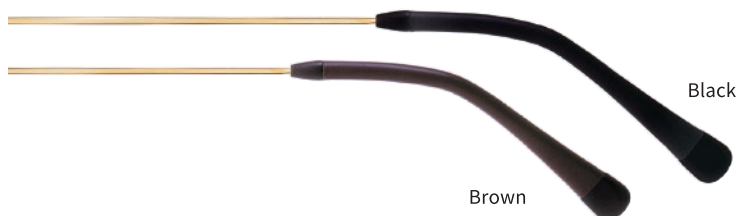

Color	Black /Brown
Hole size	1.4 ϕ — 1.5 ϕ (round hole)
Package	1) 5 pairs (Black :2 pairs / Brown :3 pairs) 2) 1 color 5 pairs (per color)
Material	Cellulose propionate, Silicone

SILICONE TEMPLE TIP (FLAT HOLE)



Color	Black / Brown / Demi
Hole size	1.4 flat — 1.5 flat
Package	1) 5 pairs (Brown:2 / Demi:2 / Black:1) 2) 1 color 5 pairs (per color)
Material	Cellulose propionate, Silicone

SILICONE TEMPLE TIP (SLIM TYPE)




Color	Black / Brown / Demi / Pearl beige / Transparent
Hole size	1.2 ϕ — 1.3 ϕ (round hole)
Package	1) 5 pairs (per color) 2) 1 color 5 pairs (per color)
Material	Cellulose propionate, Silicone



Soft material temple tip for round and flat hole

SOFURE TEMPLE TIP

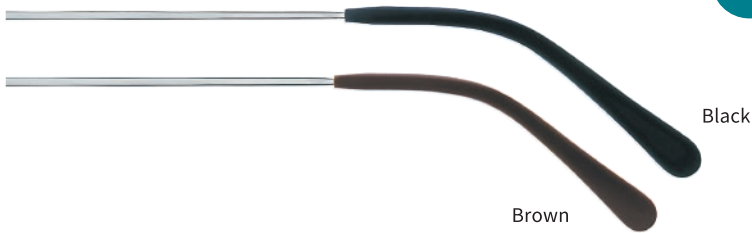


- Polyester elastomer (soft plastic).
- Very durable (cracking is also difficult to occur)
It is gentle for the skin.

Color	Black / Brown
Size	1.3φ — 1.4φ (round and flat)
Package	1) 8 pairs (Black :3 pairs/Brown :5 pairs) 2) 1 color 8 pairs (per color)
Material	Polyester elastomer


SILICONE CURL TIP FOR ADULT



SILICONE CURL TIP FOR KIDS

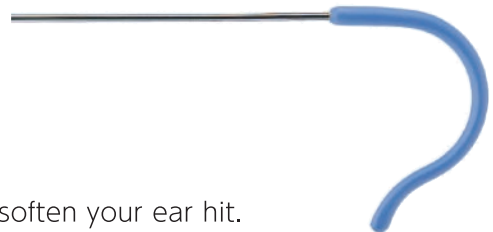

- Soft spring properties with special resin.
- Special silicone on the surface gentle to your skin and soften your ear hit.
- It can be adjusted according to the shape of the ear with hot air.
- Abundant color variation.

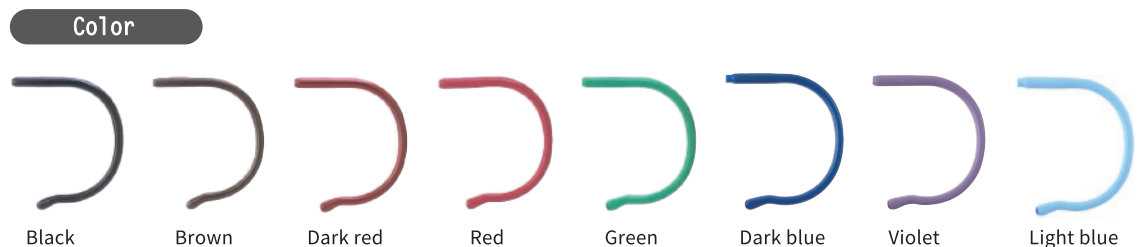

Color	Black / Brown / Dark red / Red / Green / Dark blue / Violet/ Light blue	
Size	1.3φ — 1.5φ round hole	
Package	for Adult : 1) 4 pairs (Black :1/ Brown :2 / Dark red :1)	2) 1 color 4 pairs (per color)
	for Kids : 1) 4 pairs (Light blue / Green / Red / Violet :each 1 pair)	2) 1 color 4 pairs (per color)
Material	Silicone, Polyester elastomer	



Add new function to 580 Univ. Lightest 53g ^{※1}

580 UniV[®] Pro

PD 50-80

- PD 50~80 mm (25-40 mm per right and left)
- Adjust : Temple length · Height & width & open angle of nose pad · Front angle
- With corneal vertex distance measurement scale
- 5 Lenses (Front:3pcs. & back:2pcs.)



2019

"The most emotional eyewear" Grandprix



This trial frame allows you to accurately examine progressive lenses.

< PAT no. 7132078 >

<The licence of "manufacture and sales of medical equipment" 13B3X00177>

PD	50 - 80mm (25 - 40mm per right and left)		
Nose pad	Height : 9mm (Stepless), width & open angle by 2 actions, Front and back : 5mm (Stepless)		
Temple	Length : 30mm (per 5mm) Front angle : up & down 7 degrees		
Astigmatic axis	Rotation function (5 degree)	Lens	Front : 3 pcs. / Back : 2 pc.
Weight	53g	Color	Gray

※1 Compared with the trial frame (left and right PD move independently, Nose pad moves front and back) by other companies.

※ Discontinued manufacturing : 580 UniV

Half rim type for progressive lens use

580 UniV[®] Pro-LO

- Half Rim type
- Lightest 51g ^{※1}

< PAT no. 7132078 >

<The licence of "manufacture and sales of medical equipment" 13B3X00177>

※1 Compared with the trial frame (left and right PD move independently, Nose pad moves front and back) by other companies.
 ※ Discontinued manufacturing : 580 UniV-LO






• How to handle


⚠ Turning the knob too much may cause a malfunction.

Forward and backward of nose pad

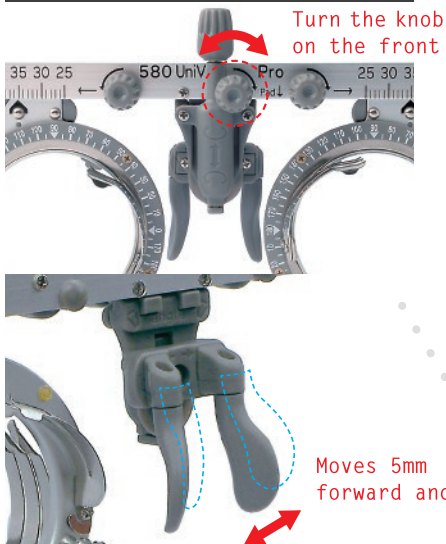


Turn the knob on the front
Nose Pad moves front and back
by 5mm.

Width & open angle of nose pad

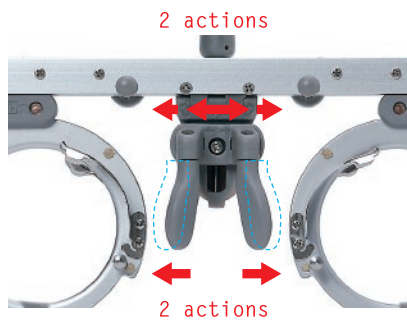


Up and down of nose pad

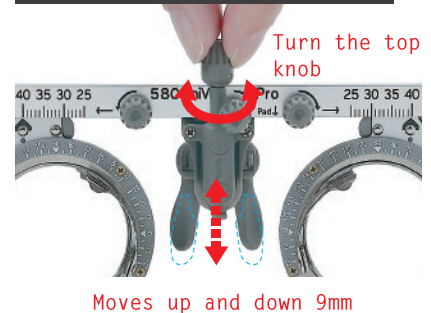


On the faces of people around the world
Versatile optometry frame that can be
matched with.

PD adjustment

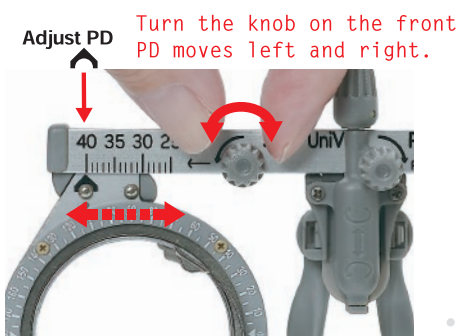

Left and right frames can be
adjusted individually.

Temple length

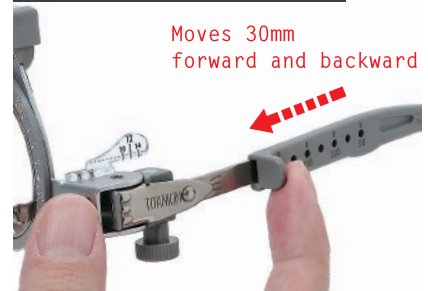


Adjust by sliding modern.

Astigmatism axis




The astigmatism axis rotates every 5 degrees.

With corneal vertex distance measurement scale

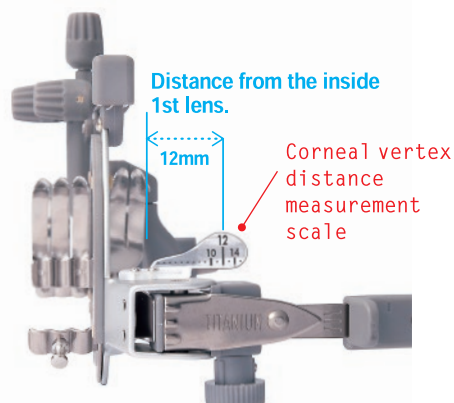

Front angle

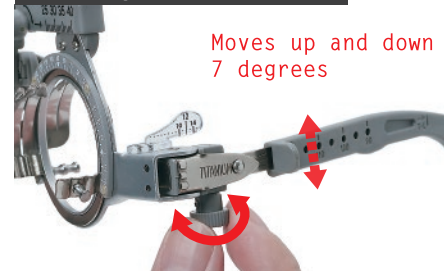


Front angle can be adjusted with a dial.



universal trial frame with one-touch PD adjustment

MC 674

PD 54 – 74

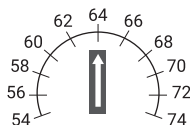

- PD 54~74 mm

<The licence of "manufacture and sales of medical equipment" 13B3X00177>

PD	54 — 74mm		
Nose pad	Hight, width, open angle by 2 actions		
Temple	Length : 30mm (per 5mm) Front angle : up & down 7 degrees		
Astigmatic axis	Rotation function (5 degree)	Lens	Front : 3 pcs. / Back : 1 pc.
Weight	50g	Color	Gray

Full rim & Straight temple

674



Half rim & Curl temple

674



Half rim & Straight temple


674



Universal trial frame for women and kids

MC 486

PD 48 – 60

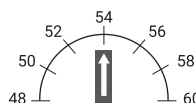

- PD 48~60 mm

<The licence of "manufacture and sales of medical equipment" 13B3X00177>

PD	48 — 60mm		
Nose pad	Hight, width, open angle by 2 actions		
Temple	Length : 20mm (per 5mm) Front angle : up & down 7 degrees		
Astigmatic axis	Rotation function (5 degree)	Lens	Front : 3 pcs. / Back : 1 pc.
Weight	48g	Color	Gray

Full rim & Straight temple

486



Half rim & Curl temple

486



Half rim & Straight temple

486



Spare for temple tip



- Straight temple tip



- Curl temple tip

Package	1 pair
For	674 series, 486 series 580 series



Standars PD fix type TRIAL FRAME

PD	54	56	58	60	62	64	66	68	70	72	74
----	----	----	----	----	----	----	----	----	----	----	----

SIMPLE BC

- Ultra lightweight (35 g) with few breakdown frames.
- Movable temple tip and gently embracing your head.
- Safe hinge structure without risk during desorption.



Simple BC for KIDS

COLOR SIMPLE BC

PD	48	50	52
----	----	----	----



PD Size	48 / 50 / 52
Temple	Length: 9mm (per 3mm) *PD52 : 12mm (per 3mm)
Astigmatic axis	Rotation function (5 degree)
Lenses	Front: 3 pcs.
Weight	33g
Color	Blue

For the professional lens use

SIMPLE BC-LO

PD	56	58	60	62	64	66	68	70	72
----	----	----	----	----	----	----	----	----	----



<The licence of "manufacture and sales of medical equipment" 13B3X00177>

PD Size	54 (simpleBC) / 56/58/60/62/64/66/68/70/72/74 (simpleBC)
Nose Pad	Height 2 action
Temple	Length 12mm (per 3mm)
Astigmatic axis	Rotation function (5 degree)
Weight	35g
Lenses	Front: 3 pcs./Back: 1 pc.
Color	Black



BLUE



BLACK

Spare temple tip


for Simple series

* Discontinued manufacturing: Pink (Color simple, spare temple tip)

other option parts for our trial frame

Clutch function for the trial frame

For examination of ptosis

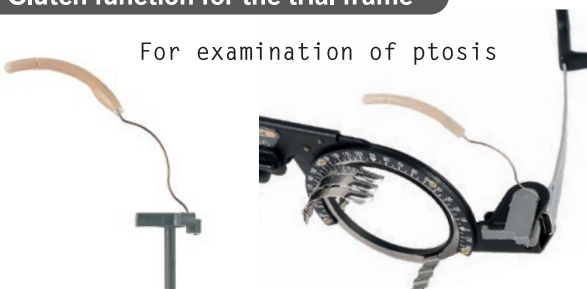

Head band

Prevents slip of the trial frame



* Some older trial frames, It may not be installed.



Tool for measuring head size

FITTING MEASURE



- For frame pre-fitting
- Temple opening width and temple tip bending length are known by mm.
- It is similar to eyeglasses.
- You can appeal your skills.


<PAT no. 3191368>

Package 1 pc.

Feel free to experience how the "high curve lens" looks

TRIAL SET FOR HIGH CURVE

- For customers considering sunglasses with high curve "Prescription glasses"

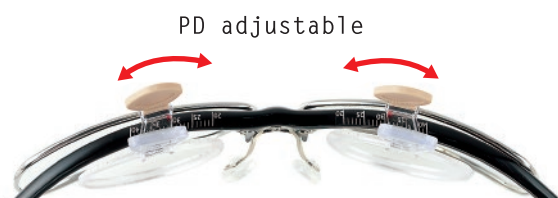

PD adjustable
Lens astigmatism inside
Spherical power lens outside

Just by attaching a high curve lens to the frame
You can experience how the high curve looks.

- It is adjustable for customer PD.
 - Put an astigmatic lens inside.
- * Since high-curve lenses are not included, please arrange at your shop.



<PAT no. 3178090>

Package

A) Frame	: 1 pc	B) Seal for lens	: 1 pair	C) Ring for astigmatic	: 1 pair
C) Ring for astigmatic	: 1 pair	D) Scale for axis degree decision	: 1 pc	E) Ring for high curve	: 3 pairs / 6 pairs / 10 pairs

For those who need help with ptosis

CLUTCH GLASS α for “EYELID PTOSIS”

alpha

Point for improvement

Durability has improved,
because by adding a one-turn twist at the base.



- Clutch glass alpha, with elastic spring, with special plastic, soft eyelid fit.
- Eliminate the discomfort of vision disorder caused by ptosis of the eyelid.
- It can be attached to a metal frame with a temple width of 5 mm or less.

⚠ Use “CLUTCH GLASS II” for plastic frame.

Package Sold separately right and left



We are accepting OEM for trial lens set.



For inquiries regarding OEM for trial lens set, please contact our sales department.



Slip off of the eyeglasses causes fatigue and eye strain.